

## Puros<sup>®</sup> Blend Particulate Allograft

Natural Blend. Proven Processing.



The Best of Both Worlds.

### 1 The Natural Mix For Bone Growth

Puros Blend Particulate Allograft is an anatomic-based mix of cortical and cancellous bone particulate which combines the space maintenance of cortical bone and the rapid remodelling of cancellous bone

### 3 A Trustworthy Source

Safety and performance you may expect from the Puros brand

Single-donor vials

Two particulate sizes to accommodate different defect sizes

### 2 Easy To Use

Retains osteoconductive properties due to the preservation of the natural bone matrix collagen and mineral composition, trabecular pattern and original porosity<sup>1</sup>

Easy handling—quick rehydration, five-year shelf life and room-temperature storage

No need to mix on site

# Puros Blend Particulate Allograft


Sterilized and preserved using the proprietary Tutoplast® process, Puros Blend Particulate Allograft is a high-quality allograft designed for large and small volume bone regeneration procedures.

## The Unique Tutoplast Process


The proprietary Tutoplast process assures the highest standard of tissue safety and quality.

The process helps preserve the valuable collagen matrix and tissue integrity while inactivating pathogens and gently removing unwanted materials, such as cells, antigens and viruses. The result is safe, biocompatible tissue.


For over 40 years, a variety of Tutoplast processed tissues have been safely used in more than five million procedures.<sup>1</sup>

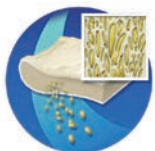


Delipidization

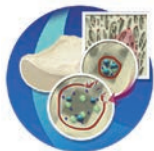


Osmotic treatment

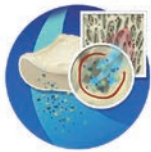


Oxidative treatment

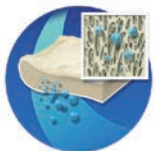

Solvent dehydration

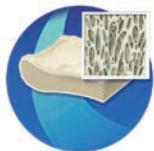

Low-dose gamma irradiation

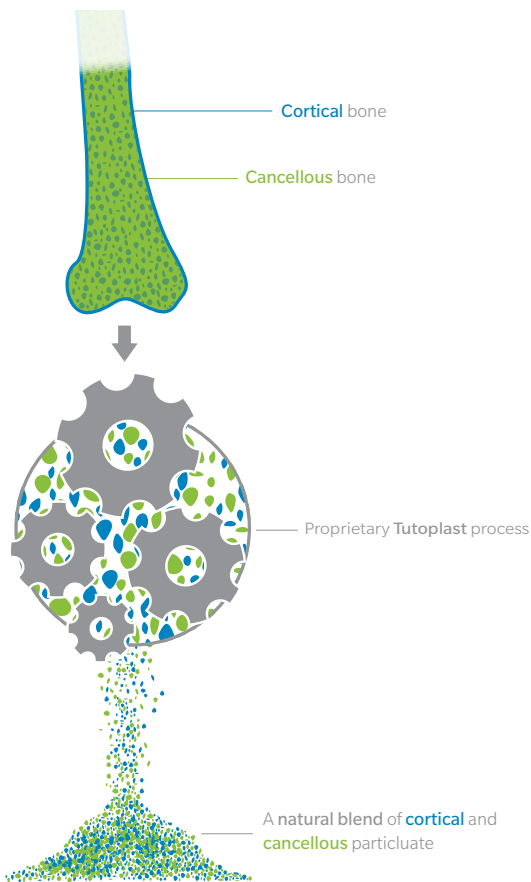

## Ordering Information

Catalog Number	Description
67800	Puros Blend Particulate, 0.5 cc, 250-1000 µm
67801	Puros Blend Particulate, 1 cc, 250-1000 µm
67802	Puros Blend Particulate, 2 cc, 250-1000 µm
67803	Puros Blend Particulate, 0.5 cc, 1000-2000 µm
67804	Puros Blend Particulate, 1 cc, 1000-2000 µm
67805	Puros Blend Particulate, 2 cc, 1000-2000 µm

Zimmer Biomet Dental offers a comprehensive line of allografts for bone augmentation needs.

<sup>1</sup>Data on file with RTI Surgical, Inc.

Contact us at 1-800-342-5454 or visit [zimmerbiometdental.com](http://zimmerbiometdental.com)

Zimmer Biomet Dental  
 Global Headquarters  
 4555 Riverside Drive  
 Palm Beach Gardens, FL 33410  
 Tel: +1-561-776-6700  
 Fax: +1-561-776-1272

All references to Zimmer Biomet Dental contained herein refer to the dental subsidiaries of Zimmer Biomet Holdings, Inc. Unless otherwise indicated, as referenced herein, all trademarks are the property of Zimmer Biomet; and all products are manufactured, distributed and marketed by Zimmer Biomet Dental (and, in the case of distribution and marketing, its authorized marketing partners). Tutoplast is a U.S. Registered trademark owned by Tutogen Medical GmbH.

Puros products are manufactured by RTI Surgical, Inc. dba Tutogen Medical GmbH. For additional product information, please refer to the individual product labeling/IFU. Product clearance and availability may be limited to certain countries/regions. This material is intended for clinicians only and does not comprise medical advice or recommendations. This material may not be copied or reprinted without the express written consent of Zimmer Biomet. ZB0396 REV A 03/18 ©2018 Zimmer Biomet. All rights reserved.

