




The single solution. Everything included.

- Multiple function Fixture Mount/Transfer minimizes chair time and financial investment.
- Efficient self-tapping flutes facilitate faster implant insertion, saving valuable surgical time.
- *Tapered SwissPlus* implant features double lead threads to provide faster insertion*.
- Narrow diameter platform with an internal hex connection is available for *Tapered SwissPlus* implant.



Preferred design and better value in a single-stage implant.

Zimmer Dental's *SwissPlus* products simplify your implant cases, while providing unmatched value. That's because each *SwissPlus* implant is packaged with a patented all-in-one Fixture Mount/Transfer which can be used for insertion, impression taking and as a customizable final abutment. Adding to the simplicity of the system, is the full assortment of *SwissPlus* prosthetics to support all possible restorative needs.

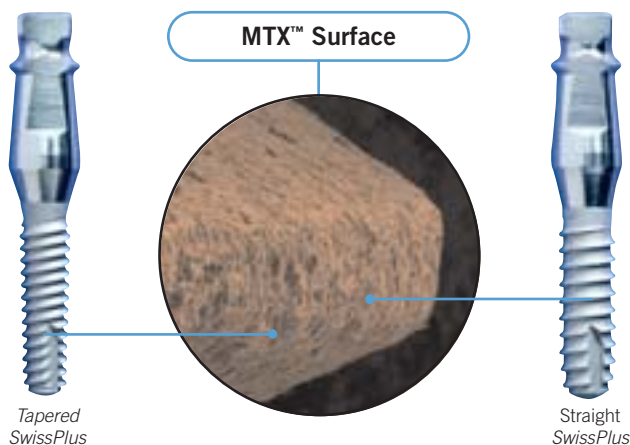

Available in our proven MTX surface**

- 1-2 micron surface roughness
- Enhanced bone apposition compared to machined titanium implants¹⁻³
- Textured with soluble, bio-compatible grit-blasting medium
- *MTX* surface processing does not round precision corners
- Cutting grooves remain intact for efficient self-tapping

* Compared to a typical single lead thread.

** Call for availability of the *MTX* surface on *SwissPlus* implants.

1. Carlsson et al, JOMI, Volume 3, number 1, pp. 21-24, 1998.
2. Buser et al, J of Biomedical Mat Res, Vol. 25, pp. 889-902, 1991.
3. Results from animal studies cannot be directly correlated to results in humans.

Multi-functional Fixture Mount/Transfer

SwissPlus implants are value-packaged with a fixture mount that functions as an impression post and preable final abutment.



Implant insertion
Fixture Mount Drill used to insert self-tapping implant.



Transfer impression
Block out the top of the transfer to make a conventional impression.



Preable final abutment
Shorten and prepare the Fixture Mount/Transfer for use as an abutment.

Prices subject to change without notice.

To learn more about this product, please visit us online at www.zimmerdental.com or to speak to a sales representative, call 800 854 7019.

To receive our eNews visit us at www.zimmerdental.com/news_eNewsLetterSignUp.aspx

For more information about our Products, Professional Programs and Continuing Education, contact us:

In the U.S. 800 854 7019
To fax an order 888 225 2483
Outside the U.S. +1 760 929 4300
Australia +61 2 4297 4182 or 1800 241 916
Canada +1 800 265 0968 or +1 905 567 2073

France +33 (0)1 45 12 35 35
Germany +49 (0)761 4584 722/723
Israel +972 (0)3 6124242
Spain +34 93 846 05 43



Zimmer Dental
1900 Aston Avenue
Carlsbad, CA 92008-7308
USA

www.zimmerdental.com