



## Zimmer® Contour Zirconia Abutments



### Strong. Beautiful. Simple.

Engineered to work with our renowned *Tapered Screw-Vent®* Implant System, the new line of Contour Zirconia Abutments from Zimmer Dental delivers an exceptional combination of strength, esthetics and simplicity. Crafted from high-strength zirconia, *Zimmer Contour Zirconia Abutments* provide a light-colored base for an esthetic, all-ceramic restoration – increasing the likelihood of overall patient satisfaction, particularly in anterior cases. These stock zirconia abutments require minimal preparation time, and with a variety of cuff heights available, there's an appropriate *Zimmer Contour Zirconia Abutment* for virtually every case.

#### Zimmer Contour Zirconia Abutments feature:

- Strong zirconia material for a natural-looking, esthetic solution
- No metal margins visible, should tissue recede
- Prepared margins to minimize preparation time
- Abutment collars that are contoured below the margin, allowing space for soft tissue around the implant
- Straight and angled designs


## A beautiful fit right off-the-shelf.

Zimmer Contour Zirconia Abutments feature a pre-defined offset margin that's lower on the buccal aspect and higher on the lingual aspect, virtually eliminating the need to prepare the abutment for a crown. The titanium interface ring provides a titanium-to-titanium connection for a precise fit between the implant and the abutment.

A recent study reported that this design may provide greater resistance to abutment micromovement and tilting during function loading, and thereby reduce the possibility of microbial leakage, screw loosening, wear and abutment fracture compared to a zirconia abutment without this seating ring<sup>1</sup>. By delivering both excellent clinical results and exceptional convenience, *Zimmer* Contour Zirconia Abutments are sure to be a beautiful fit for your patients – and your practice.



### Zimmer Contour Zirconia Abutments restoring maxillary central incisors.

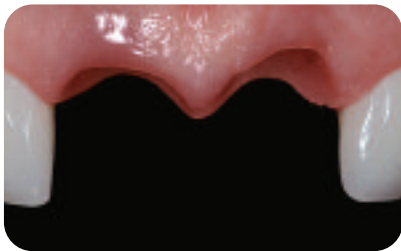

Photos ©2008 Paul S. Petrungraro, D.D.S., M.S. All rights reserved. Individual results may vary.

#### Straight



Part Number	ZRA341S	ZRA342S	ZRA451S	ZRA452S
Platform Diameter	3.5mmD	3.5mmD	4.5mmD	4.5mmD
Emergence Profile	4.5mmD	4.5mmD	5.5mmD	5.5mmD
Cuff Height	1mm	2mm	1mm	2mm

#### 17°



Part Number	ZRA351A	ZRA352A	ZRA461A	ZRA462A
Platform Diameter	3.5mmD	3.5mmD	4.5mmD	4.5mmD
Emergence Profile	5.0mmD	5.0mmD	6.0mmD	6.0mmD
Cuff Height	1mm	2mm	1mm	2mm

<sup>1</sup> Joo W. Zimmer® Contour Ceramic Abutments. Precision-engineered for strength, esthetics & clinical versatility. *Dental Inc.* 2008;2:34-38.

To learn more about *Zimmer* Contour Zirconia Abutments, please visit us online at [www.zimmerdental.com](http://www.zimmerdental.com), or to speak to a sales representative, call 1 (800) 854-7019.

To receive our eNews visit us at [http://www.zimmerdental.com/news\\_eNewsLetterSignUp.aspx](http://www.zimmerdental.com/news_eNewsLetterSignUp.aspx)

For more information about our Products and Educational Opportunities, contact us:

In the U.S. 1 (800) 854-7019  
 To fax an order 1 (888) 225-2483  
 Outside the U.S. +1 (760) 929-4300  
 Australia +61 (0)2 9950 5434 or 1 (800) 241 916  
 Canada + 1 (905) 567-2073 or 1 (800) 265-0968

France +33 (0)1 45 12 35 35  
 Germany +49 (0)761 1 56 47 0  
 Israel +972 (0)3 6124242  
 Italy +39 043 855 5573  
 Spain +34 93 846 05 43



Zimmer Dental  
 1900 Aston Avenue  
 Carlsbad, CA 92008-7308  
 USA

[www.zimmerdental.com](http://www.zimmerdental.com)