



Zimmer® Zfx™ CAD/CAM Abutments



A FITTING SOLUTION FOR YOUR PATIENTS

1 QUALITY & PRECISION

- Only source for patient-specific abutments with Zimmer's *Platform Plus™* Technology
- Comprehensive warranty on all Zfx CAD/CAM abutments
- Extensive inspection to ensure quality and fit

2 PRODUCTIVITY & PROFITABILITY

- Simplified prosthetics and streamlined workflow
- Reduced abutment inventory and associated costs
- Cost-effective alternative to traditional custom-cast abutments

3 PATIENT SPECIFIC

- Customized anatomic contours and emergence profile for improved esthetics
- Tailored restorative support for predictable performance
- Designed for proper margin placement and soft tissue management


LEVERAGE THE CLINICAL SUCCESS OF THE ZIMMER CONNECTION

Only Zimmer Zfx offers patient specific abutments with Zimmer's *Platform Plus Technology*. Incorporate the clinical benefits of Zimmer's *Platform Plus Technology* into your practice.

Zimmer Proprietary Technology – All Zimmer connections will incorporate proprietary *Platform Plus Technology* and titanium-to-titanium ring technology

Platform Plus Technology for Titanium Abutments

- Friction-Fit Connection–Virtually eliminates abutment micromovement associated with screw loosening.¹
- 1.5 mm deep Internal Hex–Shields the abutment retention screw from excess loading, distributing occlusal forces deep into the implant¹
- Lead-in Bevel–Reduces horizontal stresses better than flat “butt-joint” connections.¹⁻³

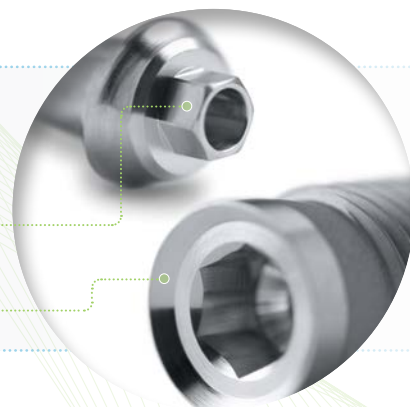

Titanium-to-Titanium Ring Technology for Zirconia Abutments

- A titanium interface ring that provides a titanium-to-titanium connection for a precise fit between the implant and abutment.⁴
- Ensures stability while distributing compressive forces more evenly throughout the implant and reducing the likelihood of component fractures.⁴

ZIMMER® ZFX™ CAD/CAM ABUTMENTS

Zimmer Connection Platform Titanium Zirconia

Zimmer Tapered Screw-Vent Implant System

	3.5mm	Titanium	Zirconia
Zimmer® Trabecular Metal™ Dental Implant	•	•	•
	4.5mm	•	•
	5.7mm	•	•

Zimmer Tapered Screw-Vent® Implant

Zimmer Screw-Vent® Implant

Note: Zirconia is available in three shades.

Visit our website for a list of competitive implant connections.

To learn more, please visit us online at www.zimmerdental.com or to speak to a sales representative, call 1 (800) 854-7019.

ORDERING AND DELIVERY

Estimated Timing* (Business Days)

	Quantity	1 unit	2-5 units	6+ units
DIGITAL FILE TRANSFER	Titanium	2-4 Days	4-6 Days	7-9 Days
	Zirconia	Call Customer Service at 1 (760) 918-3600.		
TRADITIONAL CASE (Stone Model)	Titanium	3-5 Days	5-7 Days	8-10 Days
	Zirconia	Call Customer Service at 1 (760) 918-3600.		

*Shipping time is not included. The day the order is received is not included.



Warranty

Five year comprehensive warranty on all Zfx abutments.

1. Binon PP. The evolution and evaluation of two interference-fit implant interfaces. *Postgraduate Dent.* 1996; 3:3-13.
 2. Chun HJ, Shin HS, Han CH, Lee SH. Influence of implant abutment type on stress distribution in bone under various loading conditions using finite element analysis. *Int J Oral Maxillofac Implants.* 2006; 21:105-202.
 3. Mihalko WM, May TC, Kay JF, Krause WR. Finite element analysis of interface geometry effects on the crestal bone surrounding a dental implant. *Implant Dent.* 1992;1:212-217.
 4. Won Joo, PhD, *Zimmer® Contour Ceramic Abutments: Precision-engineered for strength, esthetics and clinical versatility* ©2008 Zimmer Dental Inc.

To receive our eNews visit us at http://www.zimmerdental.com/news_eNewsLetterSignUp.aspx

For more information about our Products, Regenerative Materials and Educational Opportunities, contact us:

In the U.S. 1 (800) 854-7019
 To fax an order 1 (888) 225-2483
 Outside the U.S. +1 (760) 929-4300
 Australia +61 (0)2 9950 5434 or 1 (800) 241 916
 Canada +1 (905) 567-2073 or 1 (800) 265-0968
 Chile +562 231 5185

China + 86 21 2211 5147
 France +33 (0)1 45 12 35 35
 Germany +49 (0)761 1 56 47 0
 Israel +972 (0)3 6124242
 Italy +39 0438 37681
 Spain +34 93 846 05 43

Zimmer Dental
 1900 Aston Avenue
 Carlsbad, CA 92008-7308
 USA

www.zimmerdental.com

